

# The CapsoRetrieve Kit makes capturing the capsule simple.



1

When it's time for a bowel movement, place the Retrieval Pan on the back of the toilet bowl.



2

If needed, pour water into the pan to make it easier to see the capsule.



3

For hands-free retrieval, use the magnetic wand to capture the capsule in the pan.



4

Place the capsule into the capsule container. You can use the lid of the container to detach the capsule from the magnetic wand.



5

Close the lid and return the capsule container to your doctor's office.



OR

Or if your doctor provided you with a prepaid FedEx® envelope, place the capsule container in there and drop it off at any FedEx location.

The CapsoRetrieve Kit also contains a set of instructions that guides you through the capsule retrieval process. Read the instructions carefully and talk to your doctor if you have any questions.

## Procedure Appointment Information

Date of Procedure: \_\_\_\_\_

Day of Procedure: \_\_\_\_\_

Time of Procedure: \_\_\_\_\_

Contact Information:

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

**Caution:** Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner.

## Risk Information

Although rare, all capsule endoscopy procedures carry some risk. Possible complications include aspiration and capsule retention due to small bowel obstruction. Endoscopic placement may present additional procedure risks. Medical, endoscopic or surgical intervention may be necessary to address those complications should they occur. After ingesting the CapsoCam Plus and until the capsule is excreted, patients should not be near any strong electromagnetic fields such as MRI to avoid injury. Ask your physician for further information.



CapsoVision.com

ADV 5490 ENUS Rev E 01/23 © 2023 CapsoVision, Inc. All rights reserved.



# Capsule Endoscopy

Using CapsoCam Plus®  
The only capsule with a  
360° panoramic view


If you suffer from iron deficiency, Crohn's disease, Celiac disease, bleeding, or anemia, your doctor may want to prescribe a capsule endoscopy procedure to visualize possible abnormalities inside your small bowel (a part of your gastrointestinal (GI) tract).



Designed for Patient Freedom

## What is Capsule Endoscopy?

Capsule endoscopy is a procedure that allows your doctor to use tiny cameras inside of a capsule that is swallowed to view your small intestine as it passes through. This area is difficult to reach and see using traditional more invasive endoscopic procedures.

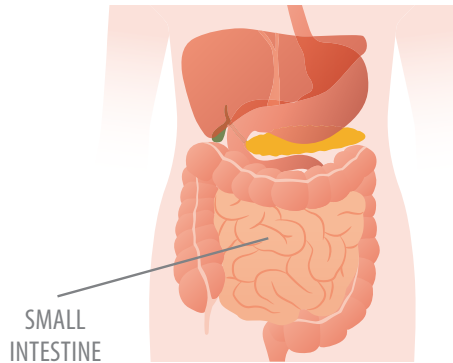

## How does CapsoCam Plus work?

CapsoCam Plus captures a unique 360° video of your digestive tract as it travels through your system. The procedure is easy:

- *Follow an easy bowel prep the day before your procedure to ensure a clean view of your bowel.*
- *Swallow the capsule with a full glass of water.*
- *Retrieve the capsule 3-30 hours later using a simple retrieval process after your bowel movement.*
- *Mail or take the capsule back to your doctor.*
- *Your doctor will review the video and share the findings with you.*



*The CapsoCam Plus capsule is approximately the size of a large multi-vitamin and can be swallowed with a few sips of water.*

## Will I feel any pain?

You should not feel any pain while swallowing the capsule or as the capsule travels through your digestive tract. The CapsoCam Plus capsule should pass naturally with a bowel movement. Unlike other endoscopic procedures, you will not require sedation.

## How do I prepare to take the capsule?

Prior to swallowing your capsule, provide your doctor with a list of all medications, allergies to medications and medical conditions. Your doctor will tell you when to start fasting. They will also provide you with pre- and post-procedure instructions. To set up the best chance for a successful outcome, carefully follow all of your doctor's instructions.

*"My patients are always surprised about how easily they can swallow the CapsoCam Plus capsule. And retrieving the capsule has been seamless for them."*

Michael Bass, M.D., Gastroenterologist  
GI Specialists of Delaware, Newark, DE

*"When given a choice, my patients prefer CapsoCam Plus. They can come in, swallow the capsule and resume their regular activities without anyone knowing that they are undergoing a procedure."*

Ian Storch, D.O., Gastroenterologist  
Gastroenterology Consultants of Long Island, New Hyde Park, NY



## How do I collect the capsule?

Before you swallow the CapsoCam Plus capsule, your doctor will provide you with the CapsoRetrieve Kit. Take the kit home and review the Instructions for Use right away. If you have any questions, be sure to talk to your doctor.

You may also want to visit the CapsoVision website where you can watch a brief video on the capsule retrieval process.



*You can also go to our website to see a short instructional video on how to retrieve the CapsoCam Plus capsule at:*

<https://capsovision.com/patient-resources/what-to-expect-with-capsocam/>

